

9853 Pacific Heights Blvd. Suite D. San Diego, CA 92121, USA Tel: 858-263-4982

Email: info@abeomics.com

36-1330: Monoclonal Antibody to Involucrin (Squamous Cell Terminal Differentiation Marker)(Clone: SY5)

Clonality: Monoclonal

Clone Name: SY5 Application: IHC Reactivity: Human Gene: IVL Gene ID: 3713 **Uniprot ID:** P07476 Format: Purified **Alternative Name:** IVI

Isotype: Mouse IgG1, kappa

Immunogen Information: Purified involucrin from human keRatinocytes

Description

It recognizes a protein of 66kDa-170kDa, identified as involucrin. In Western blotting of cultured human keratinocytes, this MAb reacts with a 120kDa protein. It stains the involucrin in a variety of sizes: 170kDa in MCF-7 cells, a doublet of ~115kDa and 150kDa in gorilla and owl monkey, 66kDa in dog, and a doublet of 105kDa in pig. Its epitope maps between codon 421-568 of human involucrin. Involucrin is expressed in a range of stratified squamous epithelia, including the cornea, which lacks a distinct cornified layer. In normal epidermis, it is first expressed in the upper spinous layers, and in keratinocyte cultures, all cells that have left the basal layer express it. Involucrin expression is altered in pathological conditions: in psoriasis and other benign epidermal hyperplasias, involucrin expression begins closer to the basal layer than normal; expression is abnormal in squamous cell carcinomas and premalignant lesions, and is reduced in severe dysplasias of the larynx and cervix.

Product Info

Amount: 100 µg

Purification: Affinity Chromatography

100 µg in 500 µl PBS containing 0.05% BSA and 0.05% sodium azide. Sodium azide is highly Content:

Store the antibody at 4°C; stable for 6 months. For long-term storage; store at -20°C. Avoid Storage condition:

repeated freeze and thaw cycles.


Application Note

Flow Cytometry (0.5-1̸µg/million cells in 0.1ml); Immunofluorescence (1-2Ã∏µg/ml); Immunohistology (Formalin-fixed) (0.1-0.2̸µg/ml for 30 min at RT); (Staining of formalin-fixed tissues is enhanced by digestion with trypsin or Protease XXV at 1mg/ml PBS for 5 min at 37 $\tilde{A}\Pi\hat{A}^{\circ}C$); (Note than enzyme digestion is better than Citrate buffer based epitope unmasking); Optimal dilution for a specific application should be determined.



9853 Pacific Heights Blvd. Suite D. San Diego, CA 92121, USA Tel: 858-263-4982

Email: info@abeomics.com

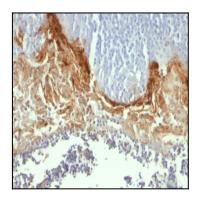

Formalin-fixed, paraffin-embedded human Tonsil stained with Involucrin Monoclonal Antibody (SY5) $\,$