

9853 Pacific Heights Blvd. Suite D. San Diego, CA 92121, USA Tel: 858-263-4982

Email: info@abeomics.com

## 36-1395: Monoclonal Antibody to Laminin, gamma 1 (LAMC1)(Clone: A5)

Clonality: Monoclonal

**Clone Name:** Α5

FACS.IF.IHC Application: Reactivity: Human, Mouse

Gene: Lamc1 **Uniprot ID:** P02468 Format: Purified **Alternative Name:** LAMC1,LAMB2 Isotype: Rat IgG2a, kappa

Immunogen Information: Murine EHS laminin prepaRation

## **Description**

Laminins are large hetero-trimeric, non-collagenous glycoproteins composed of alpha, beta, and gamma chains. This MAb reacts with laminin B2/1 chain of ~210kDa and does not cross-react with other basement membrane components or fibronectin. Its specificity was established by immunoprecipitation and immunofluorescence of human skeletal muscle and kidney with laminin chain-specific MAbs. Epithelial sheets in vivo are separated from the mesenchymal elements of the stroma by a thin layer of a specialized type of extracellular matrix termed the basement membrane (BM). This structure consists of individual components, some of which are ubiquitous in BMs and some are not. The ubiquitous ones comprise laminin (LN), entactin/nidogen (EN), collagen type IV (CIV), and large heparan sulfate proteoglycan (HSPG), which interact specifically with each other to form a continuous and regular BM. Alterations of BM integrity, from local discontinuities up to complete loss, are described in many types of human and animal epithelial neoplasms. This MAb stains uniformly all human and murine basement membranes.

## **Product Info**

Amount: 100 µg

**Purification:** Affinity Chromatography

100 µg in 500 µl PBS containing 0.05% BSA and 0.05% sodium azide. Sodium azide is highly Content:

toxic.

Store the antibody at 4°C; stable for 6 months. For long-term storage; store at -20°C. Avoid Storage condition:

repeated freeze and thaw cycles.


## **Application Note**

Flow Cytometry (0.5-1ug/million cells); Immunofluorescence (0.5-1.0ug/ml); Immunohistochemistry (Formalin-fixed) (1-2ug/ml for 30 minutes at RT)(Staining of formalin-fixed tissues requires heating tissue sections in 10mM Tris with 1mM EDTA, pH 9.0, for 45 min at 95 &degC followed by cooling at RT for 20 minutes)



9853 Pacific Heights Blvd. Suite D. San Diego, CA 92121, USA Tel: 858-263-4982

Email: info@abeomics.com

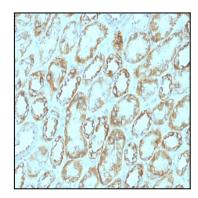

Formalin-fixed, paraffin-embedded human Renal Cell Carcinoma stained with Laminin Monoclonal Antibody (A5).