

9853 Pacific Heights Blvd. Suite D. San Diego, CA 92121, USA Tel: 858-263-4982

Email: info@abeomics.com

36-1511: Monoclonal Antibody to Neurofilament, phospho (NF-H) (Neuronal Marker)(Clone: **NE14**)

Clonality: Monoclonal **Clone Name:** NE14

Application: WB,FACS,IHC

Reactivity: Human, Mouse, Rat

Gene: NEFH Gene ID: 4744 **Uniprot ID:** P12036 Format: Purified

Alternative Name: NEFH, KIAA0845, NFH Isotype: Mouse IgG1, kappa

Immunogen Information: Crude neurofilament prepaRation from porcine spinal cord

Description

This MAb reacts with a 200kDa protein, identified as heavy sub-unit of neurofilaments (NF-H). It reacts specifically with the phosphorylated KSP/KEP segment at the C-terminus of the heavy subunit (NF-H) of neurofilaments. After dephosphorylation of neurofilaments with alkaline phosphatase, this Ab no longer binds. Neurofilaments make up the main structural elements of axons and dendrites and are found in neurons, peripheral nerves, and sympathetic ganglion cells. Neurofilaments consist of three major subunits with molecular weights of 68kDa (NF-L), 160kDa (NF-M) and 200kDa (NF-H). Anti-neurofilament stains a number of neural, neuroendocrine, and endocrine tumors. Neuromas, ganglioneuromas, gangliogliomas, ganglioneuroblastomas, and neuroblastomas stain positively for anti-neurofilament. Neurofilaments are also present in paragangliomas as well as adrenal and extra-adrenal pheochromocytomas. Carcinoids, neuroendocrine carcinomas of the skin, and oat cell carcinomas of the lung also express neurofilament.

Product Info

Amount: 100 µg

Purification: Affinity Chromatography

100 µg in 500 µl PBS containing 0.05% BSA and 0.05% sodium azide. Sodium azide is highly Content:

Store the antibody at 4°C; stable for 6 months. For long-term storage; store at -20°C. Avoid Storage condition:

repeated freeze and thaw cycles.


Application Note

Western Blot (1-2ug/ml); Flow Cytometry (1-2ug/million cells); Immunohistochemistry (Formalin-fixed) (0.25-0.5ug/ml for 30 minutes at RT)(Staining of formalin-fixed tissues requires heating tissue sections in 10mM Tris with 1mM EDTA, pH 9.0, for 45 min at 95°C followed by cooling at RT for 20 minutes);




9853 Pacific Heights Blvd. Suite D. San Diego, CA 92121, USA Tel: 858-263-4982

Email: info@abeomics.com

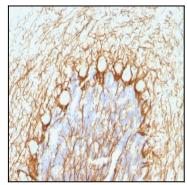

Formalin-fixed, paraffin-embedded human Cerebellum stained with Neurofilament Monoclonal Antibody (NE14).

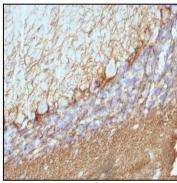

Formalin-fixed, paraffin-embedded Rat Cerebellum stained with Neurofilament Monoclonal Antibody (NE14).