

9853 Pacific Heights Blvd. Suite D. San Diego, CA 92121, USA Tel: 858-263-4982

Email: info@abeomics.com

## 36-3618: Anti-CD44v4/5 (Marker of Tumor Metastasis) Monoclonal Antibody(Clone: 3D2)

Clonality: Monoclonal

**Clone Name:** 3D2

IHC.FACS.WB Application: Reactivity: Human Gene: CD44 Gene ID: 960 **Uniprot ID:** P16070

LHR; BA-1; chondroitin sulfate proteoglycan 8 (CSPG8); Epican; Extracellular Matrix Receptor

III (ECM III); GP90 Lymphocyte Homing Adhesion Receptor; HCAM; HCELL; hematopoietic cell

**Alternative Name:** E- and L-selectin ligand; Heparan Sulfate Proteoglycan; Hermes Antigen; HAS; HUTCH I;

Hyaluronate Receptor; Indian blood group; Inlu Related p80 Glycoprotein; Ly 24; MDU2;

MDU3; MIC4; MUTCH I; Phagocytic Glycoprotein 1 (PGP-1)

Isotype: Mouse IgG1, kappa

Immunogen Information: Recombinant human full-length CD44v3-10 protein

## **Description**

This antibody recognizes an epitope encoded by exon v9 on the variant portion of human CD44. The CD44 molecule belongs to a family of cellular adhesion molecules found on a wide range of normal and malignant cells in epithelial, mesothelial and hemopoiesis tissues. CD44 is a single gene with 20 exons, of which 10 are normally expressed to encode the basic CD44 (H-CAM) molecule. The additional 10 exons (v1 to v10) are only expressed by alternative splicing of the nuclear RNA. The expression of specific cell adhesion molecule CD44 splice variants has been reported to be associated with metastasis in certain human malignancies.

## **Product Info**

Amount: 20 μg / 100 μg

200 µg/ml of Ab Purified from Bioreactor Concentrate by Protein A/G. Prepared in 10mM PBS Content:

with 0.05% BSA & 0.05% azide. Also available WITHOUT BSA & azide at 1.0mg/ml.

Antibody with azide - store at 2 to 8°C. Antibody without azide - store at -20 to -80°C. Antibody Storage condition:

is stable for 24 months. Non-hazardous.

## **Application Note**

Flow Cytometry (1-2ug/million cells); Western Blot (1-2ug/ml); Immunohistochemistry (Formalin-fixed) (1-2ug/ml for 30 minutes at RT)(Staining of formalin-fixed tissues requires heating tissue sections in 10mM Tris with 1mM EDTA, pH 9.0, for 45 min at 95°C followed by cooling at RT for 20 minutes);



9853 Pacific Heights Blvd. Suite D. San Diego, CA 92121, USA Tel: 858-263-4982

Email: info@abeomics.com

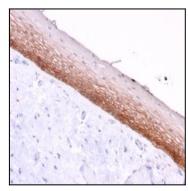


Fig. 1: Formalin-fixed, paraffin-embedded human Cervix stained with CD44v4/5 Mouse Monoclonal Antibody (3D2).

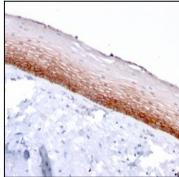


Fig. 2: Formalin-fixed, paraffin-embedded human Cervix stained with CD44v4/5 Mouse Monoclonal Antibody (3D2).

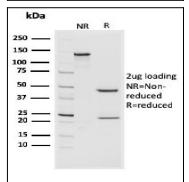


Fig. 3: SDS-PAGE Analysis Purified CD44v4/5 Mouse Monoclonal Antibody (3D2). Confirmation of Purity and Integrity of Antibody.