

9853 Pacific Heights Blvd. Suite D. San Diego, CA 92121, USA Tel: 858-263-4982

Email: info@abeomics.com

36-3793: Anti-Cytokeratin 8/18 Monoclonal Antibody(Clone: SPM192 + SPM265)

Clonality: Monoclonal

Clone Name : SPM192 + SPM265 Application : FACS,IF,WB,IHC

Reactivity: Human

 Gene :
 KRT8 & KRT18

 Gene ID :
 3856; 3875

 Uniprot ID :
 P05787; P05783

CARD2; CK8; CYK8; CYKER; Cytokeratin Endo A; DreK8; EndoA; K2C8; K8; Keratin 8; Krt 2.8;

Alternative Name: KRT8; Type-II Keratin Kb8. Cell Proliferation-inducing Gene 46 Protein; CK18;

CYK18Cytokeratin Endo B; K18; Keratin-18; Kerd; KRT18

Isotype: Mouse IgG1, kappa + Mouse IgG1, kappa

Immunogen Information: Keratin preparation from a human carcinoma (SPM192); PMC-42 human breast carcinoma

cells (SPM265)

Description

Cytokeratin 8 (CK8) belongs to the type II (or B or basic) subfamily of high molecular weight cytokeratins and exists in combination with cytokeratin 18 (CK18). This MAb cocktail recognizes all simple epithelia including glandular epithelium, for example thyroid, female breast, gastrointestinal tract, respiratory tract, and urogenital tract including transitional epithelium. All adenocarcinomas and most squamous carcinomas are positive but keratinizing squamous carcinomas are usually negative. This antibody is useful in demonstrating the presence of Paget cells; there is very little keratin 18 in the normal epidermis so only Paget cells are stained.Immunohistochemical staining with this MAb is indistinguishable from that obtained with monoclonal antibody 5D3.

Product Info

Amount: 20 μg / 100 μg

Content: 200 µg/ml of Ab Purified from Bioreactor Concentrate by Protein A/G. Prepared in 10mM PBS

with 0.05% BSA & 0.05% azide. Also available WITHOUT BSA & azide at 1.0mg/ml.

Storage condition:

Antibody with azide - store at 2 to 8°C. Antibody without azide - store at -20 to -80°C. Antibody

is stable for 24 months. Non-hazardous.

Application Note

Flow Cytometry (1-2ug/million cells); Immunofluorescence (1-2ug/ml); ,Western Blot (1-2ug/ml); ,Immunohistochemistry (Formalin-fixed) (1-2ug/ml for 30 min at RT),(Staining of formalin-fixed tissues requires heating tissue sections in 10mM Tris with 1mM EDTA, pH 9.0, for 45 min at 95°C followed by cooling at RT for 20 minutes);

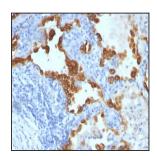


Fig. 1: Formalin-fixed, paraffin-embedded human Lung Carcinoma stained with Cytokeratin 8/18 Monoclonal Antibody (SPM192 + SPM265).