

9853 Pacific Heights Blvd. Suite D. San Diego, CA 92121, USA Tel: 858-263-4982

Email: info@abeomics.com

## 12-4087: Phospho-Chk1 (Ser317) (Clone: F10) rabbit mAb

Clonality: Monoclonal
Clone Name: Chk1S317-F10

**Application:** FACS

Reactivity: Human, Mouse
Conjugate: Unconjugated
Format: Purified

Alternative Name: Serine/threonine-protein kinase Chk1, CHK1 checkpoint homolog, Cell cycle checkpoint

kinase, Checkpoint kinase-1, CHEK1

**Isotype:** Rabbit IgG1k

Immunogen Information: A synthetic phospho-peptide corresponding to residues surrounding Ser317 of human

phospho Chk1

## **Description**

The act of DNA damaged response and cell cycle checkpoints requires the activation of four protein kinases that form the canonical ATR-Chk1 and ATM-Chk2 pathways. ATR activation requires the generation of structures containing single strand DNA (ssDNA) adjacent to double strand DNA (dsDNA). Such ssDNA is coated with replication protein A complex and attracts ATR (1,2). The accumulation of ATR to damage sites results in initial activation of ATR. ATR phosphorylates proteins at the ssDNA which are called checkpoint regulators. The accumulation and phosphorylation of these checkpoint regulators further stimulates the catalytic activity of ATR. ATR-induced Chk1 phosphorylation likely occurs at the sites of DNA damage on chromatin (3-5). The activated ATR phosphorylates Chk1 at Ser317 and Ser345 in its C-terminal regulatory domain. Phospho Chk1 is critical for DNA damage checkpoint activation, replication control, and cell viability (6-8). Functionally, ATR-mediated phosphorylation elevates phospho Chk1 catalytic activity. The N-terminal catalytic domain of Chk1 adopts an open kinase conformation and the deletion of C-terminal domain increases Chk1 catalytic activity.

## **Product Info**

**Amount :** 20 μl / 200 μl

**Content:** 1X PBS, 0.02% NaN3, 50% Glycerol, 0.1% BSA

**Storage condition:** Store at -20°C. Avoid repeated freeze and thaw cycles.

## **Application Note**

1̸µg/mL - 0.001Ã∏µg/mL. It is recommended that the reagent be titrated for optimal performance for each application. See product image legends for additional information.(0.5mg/ml)

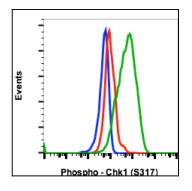


Fig-1: Flow cytometric analysis of Hela cells secondary antibody only negative control (blue) or treated imatinib (red) or treated with pervanadate (green) using Phospho-Chk1 (S317) antibody Chk1S317-F10  $0.1~\mu g/mL$ .



9853 Pacific Heights Blvd. Suite D. San Diego, CA 92121, USA Tel: 858-263-4982

Email: info@abeomics.com

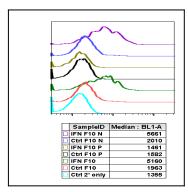


Fig 2 : Peptide blocking flow cytometric analysis of NIH3T3 cells secondary antibody only negative control (light blue) or untreated (red) or treated with IFNa + IL4 + pervanadate (green) or untreated and blocked with phospho-peptide (black) or treated and blocked with phospho peptide (gold) or untreated and blocked with non-phospho peptide (dark blue) or treated and blocked with non-phospho peptide (purple) using Phospho-Chk1 (S317) antibody Chk1S317-F10 0.1  $\mu g/mL$ .

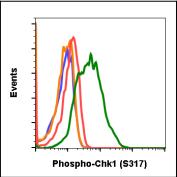


Fig-3: Flow cytometric analysis of NIH3T3 cells secondary antibody only negative control (blue) or 0.1  $\mu$ g/mL of isotype control (orange) or treated with imatinib (red) or with pervanadate (green) using Phospho-Chk1 (S317) antibody Chk1S317-F10 at 0.1  $\mu$ g/mL.