# **₩** abeomics

## 12-4347: Phospho-Akt1 (Ser473) (Clone: B9) rabbit mAb FITC conjugate

Clonality :	Monoclonal
Clone Name :	AktS473-B9
Application :	FACS
Reactivity :	Human, Mouse
Conjugate :	FITC
Format :	Conjugated
Alternative Name :	RAC-alpha serine/threonine-protein kinase, Protein kinase B, PKB
Isotype :	Rabbit IgG1k
Immunogen Information	A synthetic phospho-peptide corresponding to residues surrounding Ser473 of human phospho Akt1

### Description

Akt also known as PKB (Protein kinase B) or RAC-PK (Related to the A and C kinases) is a serine/threonine kinases that contains a pleckstrin homology (PH) domain. This protein kinase is activated by insulin and various growth and survival factors to function in a wortmannin sensitive pathway involving PI3 kinase. Akt is activated by phospholipid binding and activation loop phosphorylation at Thr308 by PDK1 and by phosphorylation within the carboxy terminus at Ser473. Phospho-Akt promotes cell survival by inhibiting apoptosis. Specifically, phospho-Akt1 has been shown to phosphorylate Bad, a member of the Bcl-2 family that promotes cell death. This phosphorylation results in the inactivation of the proapoptotic function of Bad. The Akt/phospho Akt molecule is thus considered to link extracellular survival signals (growth factors) with the apoptotic machinery (BAD). Akt is also a key mediator of the metabolic effects of insulin. Additionally, Akt has been referred to as an oncogene because it has increased activity in a number of tumors. This antibody recognizes phospho Akt phosphorylated at Ser473. This phosphorylation site is shared by all three isoforms of phospho Akt. The homologous phosphorylation sites in Akt2 and Akt3 are S474 and S472, respectively.

#### **Product Info**

Amount :	10 Tests / 100 Tests
Content :	1X PBS, 0.09% NaN3, 0.2% BSA
Storage condition :	Store at 2-8°C. Avoid repeated freeze and thaw cycles.

### **Application Note**

For flow cytometric staining, the suggested use of this reagent is 5  $\tilde{A} \square \hat{A} \mu L$  per million cells or 5  $\tilde{A} \square \hat{A} \mu L$  per 100  $\tilde{A} \square \hat{A} \mu L$  of staining volume. It is recommended that the reagent be titrated for optimal performance for each application. See product image legends for additional information.

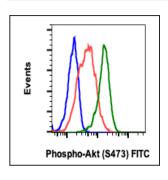


Fig-1: Flow cytometric analysis of NIH3T3 cells treated with imatinib and unstained (blue) as negative control or treated with imatinib (red) or pervadanate (PV) (green) and stained using Phospho-Akt (Ser473) antibody AKTS473-B9 FITC conjugate.

