

9853 Pacific Heights Blvd. Suite D. San Diego, CA 92121, USA Tel: 858-263-4982

Email: info@abeomics.com

30-2916: Anti-Hu CD75 PE Mab (LN1)

Clonality: Monoclonal

Clone Name : LN1
Application : FACS
Reactivity : Human
Conjugate : PE

Alternative Name: lactosamine

Isotype: Mouse IgM kappa

Immunogen Information: Stimulated human PBL

Description

Specificity: The mouse monoclonal antibody LN1 recognizes CD75, a lactosamine structure present mainly on the surface of B cell types.

CD75 is a lactosamine structure, which is present mainly on the surface of germinal center B cells. With lower level it is present on other mature B cells, and it is downregulated during differentiation to plasma cells. It is a marker which differentiates between malignant B cell types.

Product Info

Amount: 100 Tests

Purification : Purified antibody is conjugated with R-phycoerythrin (PE) under optimum conditions.

Unconjugated antibody and free fluorochrome are removed by size-exclusion chromatography.

Content: Storage Buffer: Stabilizing Tris buffered saline (TBS), pH 8.0, 15 mM sodium azide

Storage condition : Store at 2-8°C. Protect from prolonged exposure to light. Do not freeze.

Application Note

Flow cytometry: The reagent is designed for analysis of human blood cells using 10 μ l reagent / 100 μ l of whole blood or 10⁶ cells in a suspension. The content of a vial (1 ml) is sufficient for 100 tests.

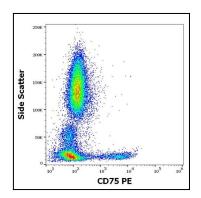


Figure 1: Flow cytometry surface staining pattern of human peripheral whole blood stained using anti-human CD75 (LN1) PE antibody (10 $\hat{l}\frac{1}{4}$ l reagent / 100 $\hat{l}\frac{1}{4}$ l of peripheral whole blood).



9853 Pacific Heights Blvd. Suite D. San Diego, CA 92121, USA Tel: 858-263-4982

Email: info@abeomics.com

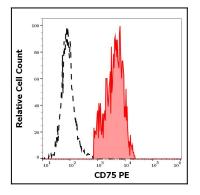


Figure 2: Separation of human CD75 positive lymphocytes (red-filled) from CD75 negative lymphocytes (black-dashed) in flow cytometry analysis (surface staining) of human peripheral whole blood stained using anti-human CD75 (LN1) PE antibody (10 \hat{l} ¹/₄l reagent / 100 \hat{l} ¹/₄l of peripheral whole blood).