

9853 Pacific Heights Blvd. Suite D. San Diego, CA 92121, USA Tel: 858-263-4982

Email: info@abeomics.com

## 30-1185PB: Pacific Blue Conjugated Anti-CD34 / Mucosialin Monoclonal Antibody (Clone:4H11[APG])

Clonality: Monoclonal **Clone Name:** 4H11[APG] Application: **FACS** Reactivity: Human Gene: CD34 Gene ID: 947 **Uniprot ID:** P28906 Mouse IgG1 Isotype:

Immunogen Information: Permanent human cell line derived from peripheral leucocytes of a patient suffering from

chronic myeloid leukaemia.

## **Description**

CD34 is a highly glycosylated monomeric 111-115 kDa surface protein, which is present on many stem cell populations. It is a well established stem cell marker, though its expression on human hematopoietic stem cells is reversible. CD34 probably serves as a surface receptor that undergoes receptor-mediated endocytosis and regulates adhesion, differentiation and proliferation of hematopoietic stem cells and other progenitors. CD34 expression is likely to represent a specific state of hematopoietic development that may have altered adhering properties with expanding and differentiating capabilities in both in vitro and in vivo conditions.

## **Product Info**

Amount: 100 tests

**Purification:** Purified antibody is conjugated with Pacific Blue NHS ester under optimum conditions and

unconjugated antibody and free fluorochrome are removed by size-exclusion chromatography.

**Content:** Formulation: Stabilizing phosphate buffered saline (PBS), pH 7.4, 15 mM sodium azide

**Storage condition :** Store at 2-8°C. Protect from prolonged exposure to light. Do not freeze.

## **Application Note**

Flow cytometry: The reagent is designed for analysis of human blood cells using 4  $\mu$ l reagent / 100  $\mu$ l of whole blood or 106 cells in a suspension. The content of a vial (0.4 ml) is sufficient for 100 tests.

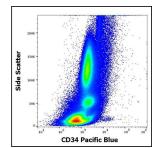


Figure 1: Flow cytometry surface staining pattern of human peripheral whole blood stained using anti-human CD34 (4H11[APG]) Pacific Blue antibody ( $4 \mu l$  reagent /  $100 \mu l$  of peripheral whole blood).



9853 Pacific Heights Blvd. Suite D. San Diego, CA 92121, USA Tel: 858-263-4982

Email: info@abeomics.com

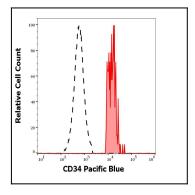


Figure 2: Separation of human CD34 positive CD45 dim stem cells (red-filled) from lymphocytes (black-dashed) in flow cytometry analysis (surface staining) of human peripheral whole blood stained using anti-human CD34 (4H11[APG]) Pacific Blue antibody (4  $\mu$ l reagent / 100  $\mu$ l of peripheral whole blood).