w abeomics

30-1185PC-Cy5.5: PerCP-Cy5.5 Conjugated Anti-CD34 / Mucosialin Monoclonal Antibody (Clone:4H11[APG])

Clonality :	Monoclonal
Clone Name :	4H11[APG]
Application :	FACS
Reactivity :	Human
Gene :	CD34
Gene ID :	947
Uniprot ID :	P28906
Isotype :	Mouse IgG1
Immunogen Information	Permanent human cell line derived from peripheral leucocytes of a patient suffering from chronic myeloid leukaemia.

Description

CD34 is a highly glycosylated monomeric 111-115 kDa surface protein, which is present on many stem cell populations. It is a well established stem cell marker, though its expression on human hematopoietic stem cells is reversible. CD34 probably serves as a surface receptor that undergoes receptor-mediated endocytosis and regulates adhesion, differentiation and proliferation of hematopoietic stem cells and other progenitors. CD34 expression is likely to represent a specific state of hematopoietic development that may have altered adhering properties with expanding and differentiating capabilities in both in vitro and in vivo conditions.

Product InfoAmount :100 testsPurified antibody is conjugated with activated tandem dye of activated Peridinin-ChlorophyllPurification :Protein-Cyanine 5.5 (PerCP-Cy5.5) under optimum conditions and unconjugated antibody and
free fluorochrome are removed by size-exclusion chromatography.Content :Formulation: Stabilizing phosphate buffered saline (PBS), pH 7.4, 15 mM sodium azide
Storage condition :Store at 2-8°C. Protect from prolonged exposure to light. Do not freeze.

Application Note

Flow cytometry: The reagent is designed for analysis of human blood cells using 4 μ l reagent / 100 μ l of whole blood or 106 cells in a suspension. The content of a vial (0.4 ml) is sufficient for 100 tests.

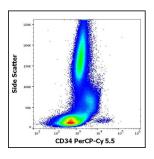


Figure 1: Flow cytometry surface staining pattern of human peripheral whole blood stained using anti-human CD34 (4H11[APG]) PerCP-Cy5.5 antibody (4 μ l reagent / 100 μ l of peripheral whole blood).

For Research Use Only. Not for use in diagnostic/therapeutics procedures.

₩ abeomics

9853 Pacific Heights Blvd. Suite D. San Diego, CA 92121, USA Tel: 858-263-4982 Email: info@abeomics.com

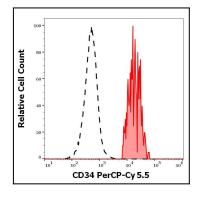


Figure 2: Separation of human CD34 positive stem cells (red-filled) from lymphocytes (black-dashed) in flow cytometry analysis (surface staining) of human peripheral whole blood stained using anti-human CD34 (4H11[APG]) PerCP-Cy5.5 antibody (4 μ l reagent / 100 μ l of peripheral whole blood).