

36-1892: Monoclonal Antibody to CD20 / MS4A1 (B-Cell Marker)(Clone : IGEL/773)

| | |
|--------------------------------|---------------------------------|
| Clonality : | Monoclonal |
| Clone Name : | IGEL/773 |
| Application : | FACS,IF,IHC |
| Reactivity : | Human |
| Gene : | MS4A1 |
| Gene ID : | 931 |
| Uniprot ID : | P11836 |
| Format : | Purified |
| Alternative Name : | MS4A1,CD20 |
| Isotype : | Mouse IgG2a, kappa |
| Immunogen Information : | Recombinant human MS4A1 protein |

Description


Recognizes a protein of 30-33kDa, which is identified as CD20. It is a non-Ig differentiation antigen of B-cells and its expression is restricted to normal and neoplastic B-cells, being absent from all other leukocytes and tissues. CD20 is expressed by pre B-cells and persists during all stages of B-cell maturation but is lost upon terminal differentiation into plasma cells. This MAAb can be used for immunophenotyping of leukemia and malignant cells, B lymphocyte detection in peripheral blood and B cell localization in tissues. It reacts with the majority of B-cells present in peripheral blood and lymphoid tissues and their derived lymphomas. In lymphoid tissue, germinal center blasts and B-immunoblasts are particularly reactive. It is a reliable antibody for ascribing a B-cell phenotype in known lymphoid tissues. Rarely, CD20-positive T-cell lymphomas have been reported. Reactivity has also been noted with Reed-Sternberg cells in cases of Hodgkin s disease, particularly of lymphocyte predominant type.

Product Info


| | |
|----------------------------|---|
| Amount : | 100 µg |
| Purification : | Affinity Chromatography |
| Content : | 100 µg in 500 µl PBS containing 0.05% BSA and 0.05% sodium azide. Sodium azide is highly toxic. |
| Storage condition : | Store the antibody at 4°C; stable for 6 months. For long-term storage; store at -20°C. Avoid repeated freeze and thaw cycles. |

Application Note

Flow Cytometry (1-2ug/million cells); Immunofluorescence (1-2ug/ml); ,Immunohistology (Formalin-fixed) (1-2ug/ml for 30 minutes at RT),(Staining of formalin-fixed tissues requires boiling tissue sections in 10mM citrate buffer, pH 6.0, for 10-20 min followed by cooling at RT for 20 minutes),

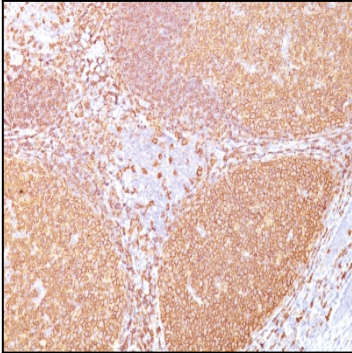

Formalin-fixed, paraffin-embedded human Tonsil stained with CD20 Monoclonal Antibody (IGEL/773)

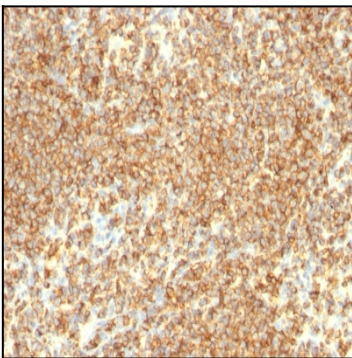

Formalin-fixed, paraffin-embedded human Lymphoma stained with CD20 Monoclonal Antibody (IGEL/773)