

9853 Pacific Heights Blvd. Suite D. San Diego, CA 92121, USA Tel: 858-263-4982

Email: info@abeomics.com

## 36-11050: Monoclonal Antibody to Plasma Cell Marker(Clone: SPM310)

Clonality: Monoclonal
Clone Name: SPM310
Application: IHC
Reactivity: Human
Format: Purified

**Isotype:** Mouse IgG2a, kappa

Immunogen Information: Pancreatic cancer-related mucin

## **Description**

It recognizes an intra-cytoplasmic antigen, which shows a very high degree of specificity for plasma cells. This antigen is present in normal as well as neoplastic plasma cells. Plasma cells, which are large lymphocytes derived from an antigen-specific B cell, secrete antibodies and are responsible for humoral immunity. Plasma cells differentiate from B cells upon stimulation by CD4+ lymphocytes. The B cell acts as an antigen-presenting cell (APC), consuming an offending pathogen, which is taken up by the B cell by phagocytosis and broken down within proteosomes. Plasma cells contain basophilic cytoplasm; their nucleus contains heterochromatin organized in a characteristic cartwheel arrangement. This MAb superbly recognizes normal and neoplastic plasma cells in routine formalin-fixed, paraffin-embedded tissue sections. It is of potential value in identifying myeloma or plasmacytoma in bone marrow or other tissues. It also helps differentiate lymphoplasmacytoid lymphoma from lymphocytic and follicular lymphoma. Note that this MAb is not suitable for staining frozen tissues.

## **Product Info**

**Amount:** 100 μg

**Purification:** Affinity Chromatography

**Content:** 100 μg in 500 μl PBS containing 0.05% BSA and 0.05% sodium azide. Sodium azide is highly


toxic.

**Storage condition :** Store the antibody at 4°C; stable for 6 months. For long-term storage; store at -20°C. Avoid

repeated freeze and thaw cycles.

## **Application Note**

Immunohistochemistry (Formalin-fixed) (0.1-0.2ug/ml for 30 minutes at RT) [(Staining of formalin-fixed tissues requires heating tissue sections in 10mM Tris with 1mM EDTA, pH 9.0, for 45 min at 95&degC followed by cooling at RT for 20 minutes);

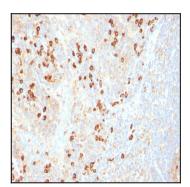

Formalin-fixed, paraffin-embedded human Tonsil stained with Plasma Cell Marker Monoclonal Antibody (SPM310).