w abeomics

14-109ACL: GATA3 Leeporter™ Luciferase Reporter-HEK293 Cell Line

Application : Functional Assay

Description

The GATA3 Leeporter[™] Luciferase Reporter cell line is a stably transfected HEK 293 cell line which expresses Renilla luciferase reporter gene under the transcriptional control of the GATA3 response element, so that the cell line is designed to measure the transcriptional activity of GATA3. As a zinc-finger transcription factor, GATA3 (GATA-binding protein 3) plays a critical role in early and late T cell differentiation, which regulates Th1/Th2 differentiation. GATA3 has been shown to induce Th2 differentiation and repress Th1 differentiation. GATA3 is also known to promote the secretion of IL-4, IL-5 and IL-13 from Th2 cells. The GATA3 induction by phorbol 12-myristate 13-acetate (PMA) is shown in Figure 1.

Product Info

Amount :	1 Vial
Content :	Each vial contains 2 ~ 3 x 10^6 cells in 1 ml of 90% FBS + 10% DMSO.
Storage condition :	Immediately upon receipt, store in liquid nitrogen.

Application Note

Application:

- Monitor the GATA3 signaling pathway activity.
- Screen for activators or inhibitors of the GATA3 signaling pathway.

Culture conditions:

Cells should be grown at 37°C with 5% CO_2 using DMEM medium (w/ L-Glutamine, 4.5g/L Glucose and Sodium Pyruvate) supplemented with 10% heat-inactivated FBS and 1% Pen/Strep, plus 3 µg/ml of Puromycin (Note: Puromycin can be omitted during the reporter cell assays).

It is recommended to quickly thaw the frozen cells upon receipt or from liquid nitrogen in a 37^oC water-bath, transfer to a tube containing 10 ml of growth medium without Puromycin, spin down cells, resuspend cells in pre-warmed growth medium without Puromycin, transfer resuspended cells to T25 flask and culture in 37^oC-CO₂ incubator.

Leave the T25 flask in the incubator for 1~3 days without disturbing or changing the medium until cells completely recover viability and become adherent. Once cells are over 90% adherent, remove growth medium and passage the cells through trypsinization and centrifugation. At first passage, switch to growth medium containing Puromycin. Cells should be split before they reach complete confluence.

To passage the cells, detach cells from culture vessel with Trypsin/EDTA, add complete growth medium and transfer to a tube, spin down cells, resuspend cells and seed appropriate aliquots of cells suspension into new culture vessels. Subcultivation ration = 1:10 to 1:20 weekly. To achieve satisfactory results, cells should not be passaged over 16 times.

For Research Use Only. Not for use in diagnostic/therapeutics procedures.



Functional validation:

A. Response of GATA3 Leeporter[™] – HEK293 cells to phorbol 12-myristate 13-acetate (PMA).

- 1. Harvest GATA3 Leeporter^M HEK293 cells and seed cells into a white solid-bottom 96-well microplate in 100 μ l of growth medium at 5 x 10⁴ cells/well.
- 2. Incubate cells at 37° C in a CO₂ incubator for overnight.
- 3. The next day, stimulate cells with different concentrations of PMA.
- 4. Incubate at 37° C in a CO₂ incubator for 16 hours.
- 5. Equilibrate the plate to room temperature for 10 minutes.

6. Add 50 μ l of luciferase assay reagent (Abeomics, Cat #17-1101; Refer to the reagent datasheet for the detailed luciferase assay protocol) per well.

7. Read the plate in 1-5 minutes to measure luminescence using a microplate luminometer.

LIMITED USE RESTRICTIONS:

THIS PRODUCT IS SOLELY FOR IN VITRO RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR THERAPEUTIC USE.

By use of this product, user agrees to be bound by the terms of this limited use statement.

This product is <u>solely for Internal Research Purposes</u> and <u>not for Commercial Purposes</u>. Commercial Purposes include, but are not limited to (1) use of the cell line in manufacturing; (2) use of the cell line to provide a service, information or data; (3) use of the cell line for therapeutic, diagnostic or prophylactic purposes; or (4) resale of the cell line whether or not such cell lines are resold for use in research. <u>The buyer cannot sell, give or otherwise transfer this product to a third party.</u>

Commercial License Agreement is available for non-research use if applicable. Please contact Abeomics (info@abeomics.com).

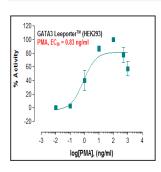


Fig-1:. Induction of GATA3 activity by PMA in GATA3 Leeporter[™] - HEK293 cells.

