

9853 Pacific Heights Blvd. Suite D. San Diego, CA 92121, USA Tel: 858-263-4982

Email: info@abeomics.com

36-1029: Monoclonal Antibody to Carcinoembryonic Antigen (CEA) / CD66(COL-1 + CEA31 + C66/261)

Clonality: Monoclonal

Clone Name: COL-1 + CEA31 + C66/261

Application: FACS, IF, IHC Reactivity: Human Gene: CEACAM5 Gene ID: 1048 **Uniprot ID:** P06731 Format: Purified **Alternative Name:** CEACAM5.CEA

Isotype: Mouse/IgG2a (COL-1) + Mouse/IgG1 (CEA31 and C66/261)

Human colon carcinoma extract (COL-1 & CEA31); Recombinant full-length human CEA protein Immunogen Information:

(C66/261)

Description

This antibody recognizes proteins of 80-200kDa, identified as different members of CEA family. CEA is synthesized during development in the fetal gut and is re-expressed in increased amounts in intestinal carcinomas and several other tumors. This MAb does not react with nonspecific cross-reacting antigen (NCA) and with human polymorphonuclear leucocytes. It shows no reaction with a variety of normal tissues and is suitable for staining of formalin/paraffin tissues. CEA is not found in benign glands, stroma, or malignant prostatic cells. Antibody to CEA is useful in detecting early foci of gastric carcinoma and in distinguishing pulmonary adenocarcinomas (60-70% are CEA+) from pleural mesotheliomas (rarely or weakly CEA+). Anti-CEA positivity is seen in adenocarcinomas from the lung, colon, stomach, esophagus, pancreas, gallbadder, urachus, salivary gland, ovary, and endocervix.

Product Info

Amount: 100 µg

Purification: Affinity Chromatography

100 µg in 500 µl PBS containing 0.05% BSA and 0.05% sodium azide. Sodium azide is highly Content:

Store the antibody at 4°C; stable for 6 months. For long-term storage; store at -20°C. Avoid Storage condition:

repeated freeze and thaw cycles.


Application Note

Flow Cytometry (0.5-1ug/million cells); Immunofluorescence (0.5-1ug/ml); ,Immunohistochemistry (Formalin-fixed) (1-2ug/ml for 30 minutes at RT),(Staining of formalin-fixed tissues requires heating tissue sections in 10mM Tris with 1mM EDTA, pH 9.0, for 45 min at 95 °C followed by cooling at RT for 20 minutes),



9853 Pacific Heights Blvd. Suite D. San Diego, CA 92121, USA Tel: 858-263-4982

Email: info@abeomics.com

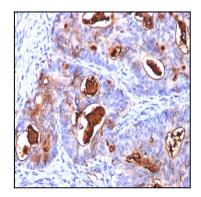

Formalin-fixed, paraffin-embedded human Colon Carcinoma stained with CEA Monoclonal Antibody (COL-1+CEA31+C66/261)