

36-1197: Monoclonal Antibody to Blood Group Antigen B (CD173)(Clone : HEB-20)

Clonality :	Monoclonal
Clone Name :	HEB-20
Application :	IF,IHC
Reactivity :	Human
Gene :	ABO
Gene ID :	28
Uniprot ID :	P16442
Format :	Purified
Alternative Name :	ABO
Isotype :	Mouse IgG1, kappa
Immunogen Information :	Mixture of erythrocytes of group B and glycoprotein fraction isolated from saliva of secretors with blood group B

Description


The antibody HEB-20 reacts with human blood group B. The specificity of the antibody HEB-20 was confirmed by comparison of specificity and reactivity to standard reagent using >5.000 samples of blood. The MAb HEB-20 shows specific staining of erythrocytes and vascular epithelium of blood group B controls and no staining in group A controls. This MAb is applicable for tissue staining in tumor patients with blood groups B and AB. Blood group antigens are generally defined as molecules formed by sequential addition of saccharides to the carbohydrate side chains of lipids and proteins detected on erythrocytes and certain epithelial cells. The A, B and H antigens are reported to undergo modulation during malignant cellular transformation. Blood group related antigens represent a group of carbohydrate determinants carried on both glycolipids and glycoproteins. They are usually mucin type, and are detected on erythrocytes, certain epithelial cells, and in secretions of certain individuals. Sixteen genetically and biosynthetically distinct but inter related specificities belong to this group of antigens, including A, B, H, Lewis A, Lewis B, Lewis X, Lewis Y, and precursor type 1 chain antigens.

Product Info

Amount :	100 µg
Purification :	Affinity Chromatography
Content :	100 µg in 500 µl PBS containing 0.05% BSA and 0.05% sodium azide. Sodium azide is highly toxic.
Storage condition :	Store the antibody at 4°C; stable for 6 months. For long-term storage; store at -20°C. Avoid repeated freeze and thaw cycles.

Application Note

Agglutination; ,Immunofluorescence (2-4ug/ml); ,Immunohistology (Formalin-fixed) (1-2ug/ml for 30 minutes at RT),(Staining of formalin-fixed tissues requires boiling tissue sections in 10mM citrate buffer, pH 6.0, for 10-20 min followed by cooling at RT for 20 minutes),

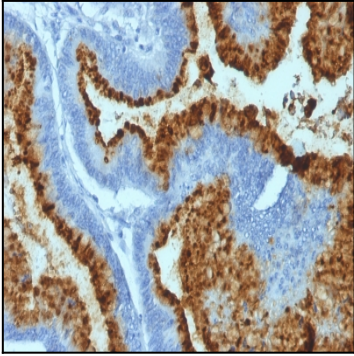

Formalin-fixed, paraffin-embedded human Colon Carcinoma stained with Blood Group B Monoclonal Antibody (HEB-20).